

# Amendia

The Elements of Healing

## ZEUS-T

TRANSFORAMINAL LUMBAR INTERBODY

Z<sub>A</sub>

Z<sub>C</sub>

Z<sub>L</sub>

Z<sub>O</sub>

Px

L

Sp

Sa

Dd

Cm

Amendia

Z<sub>T</sub>

ZEUS TRANSFORAMINAL  
LUMBAR INTERBODY



Sv<sub>T</sub>

Z<sub>P</sub>

Z<sub>V</sub>

## SURGICAL TECHNIQUE GUIDE

# Zeus-T

## Transforaminal Lumbar Interbody Device

### Table of Contents

Features and Benefits	1
Instrument Guide	2-5
Implant Guide	6-7
Surgical Technique Guide	8-14
Important Product Information	15-18



**1755 West Oak Parkway**  
**Marietta, GA 30062**  
**Phone: 877-755-3329**  
**Fax: 877-420-1213**  
**[www.Amendia.com](http://www.Amendia.com)**  
**[info@amendia.com](mailto:info@amendia.com)**

**Disclaimer:**

The surgical technique shown is for illustrative purposes only. The technique(s) actually employed in each case will always depend upon the medical judgment of the surgeon before and during surgery as to the best mode of treatment for each patient. Please reference the 510K or package insert for additional information and a complete list of intended indications, warnings, precautions, and other medical information.

## Features & Benefits

- Constructed of PEEK
- Tantalum radiopaque markers to optimize visibility and placement
- Angled teeth prevent implant migration
- Large graft cavity to facilitate fusion
- Design ensures a high degree of compressive strength and dimensional stability
- Bullet Nose option has tapered leading edge that distracts endplates and facilitates delivery



Available Sizes (mm)	
10 W x 28 L x 7 H	11 W x 36 L x 7 H
10 W x 28 L x 8 H	11 W x 36 L x 8 H
10 W x 28 L x 9 H	11 W x 36 L x 9 H
10 W x 28 L x 10 H	11 W x 36 L x 10 H
10 W x 28 L x 11 H	11 W x 36 L x 11 H
10 W x 28 L x 12 H	11 W x 36 L x 12 H
10 W x 28 L x 13 H	11 W x 36 L x 13 H
10 W x 28 L x 14 H	11 W x 36 L x 14 H
10 W x 28 L x 15 H	11 W x 36 L x 15 H



Available Sizes (mm)
11 W x 28 L x 7 H
11 W x 28 L x 8 H
11 W x 28 L x 9 H
11 W x 28 L x 10 H
11 W x 28 L x 11 H
11 W x 28 L x 12 H
11 W x 28 L x 13 H
11 W x 28 L x 14 H
11 W x 28 L x 15 H
11 W x 36 L x 7 H
11 W x 36 L x 8 H
11 W x 36 L x 9 H
11 W x 36 L x 10 H
11 W x 36 L x 11 H
11 W x 36 L x 12 H
11 W x 36 L x 13 H
11 W x 36 L x 14 H
11 W x 36 L x 15 H

# Instrument Guide




Inserter  
LM-10075



Curved Tamp  
VLM-10022



Double-Handle Curved Rasp  
VLM-10023

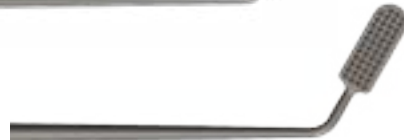

Curved Chisel  
VLM-10053



# Instrument Guide



T-Handle, Hudson Quick Connect  
LM-10019



Double-Handle Curved Oblong Curette  
VLM-10024



Double-Handle Curved Round Curette  
VLM-10025




Double-Handle Curved Box Curette  
VLM-10026

# Instrument Guide



Curved Trial

Part No.	Size
LM-10005	7 mm
LM-10080	8 mm
LM-10006	9 mm
LM-10081	10 mm
LM-10007	11 mm
LM-10017	12 mm
LM-10008	13 mm
LM-10083	14 mm
LM-10015	15 mm
LM-10082	16 mm

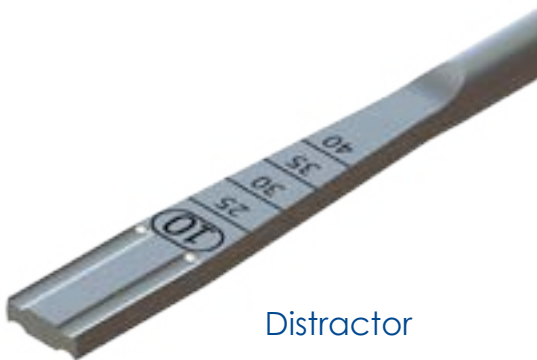

Distractor

Part No.	Size
LM-10049	7 mm
LM-10070	8 mm
LM-10050	9 mm
LM-10071	10 mm
LM-10051	11 mm
LM-10072	12 mm
LM-10052	13 mm
LM-10073	14 mm
LM-10078	15 mm
LM-10074	16 mm

# Instrument Guide

## Bullet Nose



Curved Trial

Part No.	Size
67-281107T	7 mm
67-281108T	8 mm
67-281109T	9 mm
67-281110T	10 mm
67-281111T	11 mm
67-281112T	12 mm
67-281113T	13 mm
67-281114T	14 mm
67-281115T	15 mm

Part No.	Size
67-361107T	7 mm
67-361108T	8 mm
67-361109T	9 mm
67-361110T	10 mm
67-361111T	11 mm
67-361112T	12 mm
67-361113T	13 mm
67-361114T	14 mm
67-361115T	15 mm



Ratcheting Inserter  
67-01  
(Available Upon Request)

# Implant Guide

## Available Sizes

Part No.	Size
1001107C	11 W x 28 L x 7 H
1001108C	11 W x 28 L x 8 H
1001109C	11 W x 28 L x 9 H
1001110C	11 W x 28 L x 10 H
1001111C	11 W x 28 L x 11 H
1001112C	11 W x 28 L x 12 H
1001113C	11 W x 28 L x 13 H
1001114C	11 W x 28 L x 14 H
1001115C	11 W x 28 L x 15 H
100C1107	11 W x 36 L x 7 H
100C1108	11 W x 36 L x 8 H
100C1109	11 W x 36 L x 9 H
100C1110	11 W x 36 L x 10 H
100C1111	11 W x 36 L x 11 H
100C1112	11 W x 36 L x 12 H
100C1113	11 W x 36 L x 13 H
100C1114	11 W x 36 L x 14 H
100C1115	11 W x 36 L x 15 H





# Implant Guide

## Available Sizes

### Bullet Nose

Part No.	Size
67-281107	10 W x 28 L x 7 H
67-281108	10 W x 28 L x 8 H
67-281109	10 W x 28 L x 9 H
67-281110	10 W x 28 L x 10 H
67-281111	10 W x 28 L x 11 H
67-281112	10 W x 28 L x 12 H
67-281113	10 W x 28 L x 13 H
67-281114	10 W x 28 L x 14 H
67-281115	10 W x 28 L x 15 H
67-361107	11 W x 36 L x 7 H
67-361108	11 W x 36 L x 8 H
67-361109	11 W x 36 L x 9 H
67-361110	11 W x 36 L x 10 H
67-361111	11 W x 36 L x 11 H
67-361112	11 W x 36 L x 12 H
67-361113	11 W x 36 L x 13 H
67-361114	11 W x 36 L x 14 H
67-361115	11 W x 36 L x 15 H



# Surgical Technique Guide

## Step 1: Make Incision

Incise and identify anatomical landmarks. Locate the facets, pars interarticularis, lamina, spinous processes, and transverse processes (**Figure 1**).



Figure 1

## Step 2: Distraction

Place a lamina spreader at the base of the spinous processes of the appropriate levels and apply distraction. This maneuver temporarily opens the posterior disc space and promotes increased exposure for both decompression and delivery of the **Implant**. Alternatively, pedicle screws may be implanted and used for distraction (**Figure 2**).



Figure 2

## Step 3: Create the Transforaminal Window

Prepare a window for a transforaminal approach using the **Curved Chisel** to remove the inferior articular process of the cranial vertebra and the superior articular process of the caudal vertebra (**Figure 3**).



Figure 3

## Step 4: Discectomy

Remove disc material from the intervertebral disc space using the **Curved Curettes**. The use of **Curved Curettes** will help facilitate removal of material in the far lateral disc space (**Figure 4**). Care must be taken as the anterior and lateral walls of the annulus must be preserved to provide additional support for the **Cage**. If necessary, additional distraction may also be applied at this time.

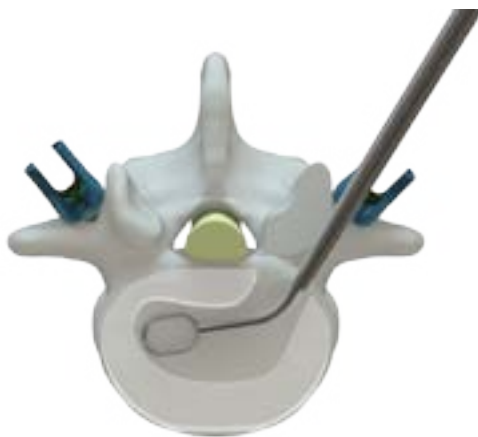


Figure 4

Remove the cartilaginous layers from the surface of the vertebral endplates with a **Double Handled Curved Rasp** until bleeding bone is attained. Sufficient cleaning of the endplates is essential for vascular supply to the bone graft; yet excessive cleaning could damage and weaken the endplate.

# Surgical Technique Guide

## STEP 5: Trial Selection

Connect an appropriately sized **Zeus-T Curved Cage Trial** to the **T-Handle**. Carefully insert the selected **Trial** via the transforaminal window into the disc space (**Figure 5**), applying gentle impaction. Fluoroscopy can assist in confirming the fit and geometry of the **Trial**.

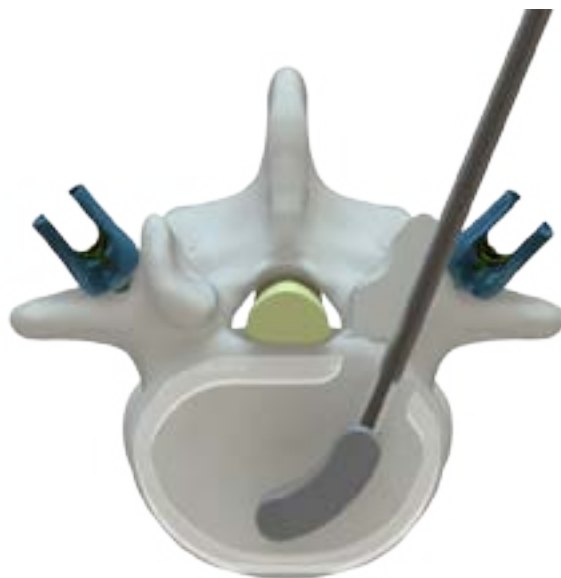


Figure 5

With the segment fully distracted, the **Trial** must fit tightly and accurately between the endplates in order to ensure that disc height will be preserved when the distraction is released.

If the **Trial** appears too small or too tight, try the next larger or smaller size until the appropriate fit is achieved. Using the largest possible **Implant** maximizes segment stability by creating tension on the longitudinal ligament and the annulus fibrosus. When the correct **Zeus-T Curved Cage** size has been determined, distraction can be temporarily released to verify proper fit.

*Note: The **Trials** are not for implantation and must be removed prior to insertion of the **Cage**.*

# Surgical Technique Guide

## STEP 6: Implant Insertion

Select the **Zeus-T Curved Cage** corresponding to the **Trial** size determined in step 5. Attach the **Cage** onto the **Implant Inserter** and fill the **Cage** with autologous bone. When the **Cage** is ready for implantation, distract the segment again.

*Note: After trialing but before the **Zeus-T Cage** is implanted, the anterior and lateral disc space should be filled with autologous bone.*

Ensure that the orientation of the **Implant** is correct and then insert the **Cage** into the intervertebral disc space. Slight impaction on the **Implant Inserter** may be necessary (**Figure 6**).

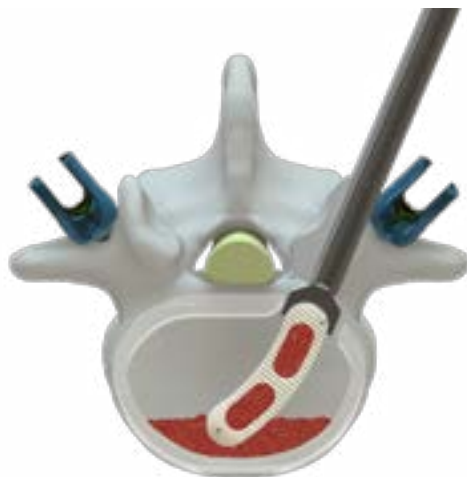


Figure 6

# Surgical Technique Guide

The optimal placement of the **Cage** is in the anterior half of the intervertebral space. Depending on the size of the vertebrae, the anterior rim of the **Cage** will be 6–8 mm posterior to the anterior edges of adjacent vertebral bodies (**Figure 7**). To help facilitate **Implant** placement, either the **Curved Tamp** can be used, or the center shaft of the **Insertor** can be removed and used in the same fashion. Use radiographic imaging to verify the A/P position of the **Cage** relative to the vertebral bodies. Once in position, remove all instruments and fill the remaining disc space with bone graft (**Figure 8**).



Figure 7



Figure 8

\*Use radiographic imaging to verify the A/P and Lateral position of the **Bullet Nose Cage** relative to the vertebral bodies (**Figure 9-A, 9-B**).

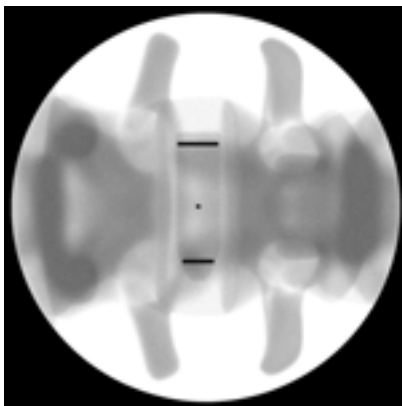


Figure 9-A

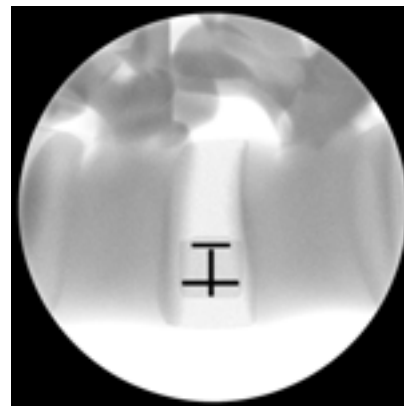


Figure 9-B

## STEP 7: Supplemental Fixation

Because additional posterior fixation considerably enhances the biomechanical stability of the motion segment as well as the stability of the device, the **Zeus-T Cage** must be used in conjunction with supplemental posterior fixation.

Retract the muscle layer and insert pedicle screws of choice if not previously implanted. The final steps of the fixation procedure (e.g. rod insertion, tightening, compression) are completed after implantation of the **Cage**. Please refer to the specific system's surgical technique manual for user instructions.

## STEP 8: Revision Surgery (implant removal)

If removal of the **Implant** is desired at a later time, then the entire posterior surface of the **Implant** should be exposed. The **Insertor** is threaded securely into the **Implant** and is removed by pulling the **Insertor** in an upward motion away from the spine.



# Important Product Information

CAUTION: Federal law (USA) restricts these devices to sale by or on the order of a physician. Implants and disposable instruments single use only.

## Description:

The Zeus® Intervertebral Body Fusion Devices are used to maintain disc space distraction in skeletally mature adults requiring intervertebral body fusion. They are designed to be used in conjunction with supplemental spinal fixation instrumentation. The series is comprised of cages of various fixed heights and shapes for placement in the cervical or lumbar spine. There are different cages designed for specific regions of the spine and approaches to the spine. Each cage has a hollow center to allow placement of autograft inside of the cage. Ridges on the superior and inferior surfaces of the device help to grip the endplates and prevent expulsion.

## Indications for Use:

Zeus® Cervical (Zeus-C®) Intervertebral Body Fusion Devices are indicated for use in skeletally mature patients with degenerative disc disease (DDD) of the cervical spine with accompanying radicular symptoms at one disc level. DDD is defined as discogenic pain with degeneration of the disc confirmed by patient history and radiographic studies. Zeus-C® cages are used to facilitate intervertebral body fusion in the cervical spine and are placed via an anterior approach at the C3 to C7 disc levels using autograft bone. Zeus-C® cages are to be used with supplemental fixation. Patients should have at least six (6) weeks of non-operative treatment prior to treatment with an intervertebral cage.

The Zeus® Lumbar Intervertebral Body Fusion Devices are indicated for intervertebral body spinal fusion procedures in skeletally mature patients with degenerative disc disease (DDD) at one level from L2-S1. DDD is defined as discogenic back pain with degeneration of the disc confirmed by patient history and radiographic studies. These DDD patients may also have up to Grade I spondylolisthesis or retrolisthesis at the involved level. Zeus® implants are to be used with autogenous bone graft and implanted via an anterior (Zeus-A®), transforaminal (Zeus-T®), versatile (Zeus-V®), Zeus-L®, Zeus-O®) or posterior (Zeus-P®) approach. The Zeus® lumbar cages are to be used with supplemental fixation. Patients should have at least six (6) months of non-operative treatment prior to treatment with an intervertebral cage.

## Warning:

An entirely satisfactory result is not always achieved in every surgical case. This particularly applies to spinal surgery, in which numerous external factors may compromise the results. Patients with previous spinal surgery at the levels to be treated may have different clinical outcomes compared to those without a previous surgery. The risk of device expulsion and migration is higher without the use of supplemental fixation.

## Contraindications include, but are not limited to:

- Spondylolisthesis higher than grade I
- Reduced bone density, which does not guarantee a sufficient resting stability (e. g. osteoporosis)
- Fractures
- Tumors
- Scoliosis
- Active infection
- Allergy to tantalum or PEEK
- Signs of local inflammation
- Fever or leukocytosis
- Morbid obesity
- Pregnancy
- Mental illness
- Suspected or documented allergy or intolerance to composite materials
- Any case not needing a fusion
- Any case not described in the indications
- Any patient unwilling to cooperate with postoperative instructions
- Patients with a known hereditary or acquired bone friability or calcification problem should not be considered for this type of surgery
- These devices must not be used for pediatric cases, nor where the patient still has general skeletal growth
- Spondylolisthesis unable to be reduced to Grade 1
- Any case where the implant components selected for use would be too large or too small to achieve a successful result
- Any case that requires the mixing of metals from two different components or systems
- Any patient having inadequate tissue coverage over the operative site or inadequate bone stock or quality
- Any patient in which implant utilization would interfere with anatomical structures or expected physiological performance
- Prior fusion at the level to be treated
- Any other condition which would preclude the potential benefit of spinal implant surgery, such as the presence of tumors or congenital abnormalities, fracture local to the operating site, elevation of sedimentation rate unexplained by other diseases, elevation of white blood count (WBC), or a marked left shift in the WBC differential count

## Additional contraindications for the Zeus-L® IBF Device include but are not limited to:

- Symptomatic level at L5-S1
- Lumbar deformities with more than 30° of rotation
- Retroperitoneal scarring on both left and right sides (e.g. due to abscess or prior surgery)
- Need for direct nerve decompression through the same approach.

# Important Product Information

## Precautions:

Intervertebral body fusion should only be undertaken after the surgeon has had hands-on training in these methods of spinal fixation, and has become thoroughly knowledgeable about spinal anatomy and biomechanics. Surgical technique manuals are available for detailed instructions on the correct use of the Zeus® cages. The contents of these manuals alone are not adequate for complete instruction in the use of this system. Even for surgeons already experienced in spinal instrumentation and intervertebral body fusion procedures, new skills may be required that are best learned by working with an experienced surgeon or through a course of formal instruction with laboratory training. Lack of experience or expertise with these implants may result in complications. Patients with previous spinal surgery at the level(s) to be treated may have different clinical outcomes compared to those without a previous surgery.

## Possible Adverse Effects:

Possible adverse events or complications associated with the Zeus® Intervertebral Body Fusion devices may include, but are not limited to:

- Bending or fracture of implant. Loosening of the implant.
- Implant material sensitivity, or allergic reaction to a foreign body.
- Infection, early or late.
- Decrease in bone density due to stress shielding.
- Pain, discomfort, or abnormal sensations due to the presence of the device.
- Nerve damage due to surgical trauma or presence of the device. Neurological difficulties including bowel and/or bladder dysfunction, impotence, retrograde ejaculation, radicular pain, tethering of nerves in scar tissue, muscle weakness, and paraesthesia.
- Vascular damage could result in catastrophic or fatal bleeding. Malpositioned implants adjacent to large arteries or veins could cause erosion of these vessels and catastrophic bleeding in the later postoperative period.
- Dural tears experienced during surgery could result in need for further surgery for dural repair, a chronic CSF leak or fistula, and possible meningitis.
- Bursitis.
- Paralysis.
- Death.
- Spinal cord impingement or damage.
- Fracture of bony structures.
- Reflex sympathetic dystrophy.
- If a pseudarthrosis occurs coupled with the Zeus® cages, a mechanical grinding action could possibly occur which might generate wear debris. Most types of wear debris have shown the potential of initiating local osteolysis in articulating joints.
- Degenerative changes or instability in segments adjacent to fused vertebral levels.

Additional surgery may be necessary to correct some of these adverse events.

## Material Specification:

The Zeus® Intervertebral Body Fusion Devices are made from the PEEK (polyetheretherketone) radiolucent material with embedded tantalum x-ray markers as specified in ASTM F2026 and ASTM F560, respectively. No warranties, expressed or implied, are made.

## Packaging:

Packages for each of the components should be intact upon receipt. Damaged packages and products should not be used and should be returned to AMENDIA.

## Sterilization:

Products not clearly marked as sterile should be assumed non-sterile.

### For Sterile Implants and Instruments:

Implants and instruments provided sterile will be clearly labeled as such in an unopened sterile package provided by AMENDIA. The contents are sterile unless the package is damaged, opened, or the expiration date on the device label has passed. The integrity of the packaging should be checked to ensure that the sterility of the contents is not compromised. Implants supplied sterilized from AMENDIA must not be re-sterilized.

### For Non-Sterile Implants and Instruments:

Implants and instruments used in surgery not clearly labeled as sterile must be sterilized by the hospital prior to use. Remove all packaging materials prior to sterilization where applicable.

Only sterile products should be placed in the operative field.

## Product Complaints:

Any health care professional (e.g. customer or user) who has experienced dissatisfaction in the services of AMENDIA or who has any complaints about AMENDIA products referring to quality, identity, durability, reliability, safety, effectiveness, and/or performance, should notify this to the sales representative, distributor, or AMENDIA customer service. Further, if any of the devices, instruments or components ever malfunction, (i.e. do not meet any of their performance specifications or otherwise do not perform as intended), or are suspected of doing so, the distributor should be notified immediately. If any AMENDIA product ever "malfunctions" and may have caused or contributed to the death or serious injury of a patient, the distributor should be notified immediately by telephone, FAX or written correspondence. When filing a complaint, please provide the component(s) name and number, lot number(s), your name and address, the nature of the complaint and notification of whether a written report from the distributor is requested.

# Important Product Information

Manufacturer: AMENDIA, 1755 West Oak Parkway, Marietta, GA 30062, 877-755-3329 (Toll Free), 770-575-5200 (Main), 877-420-1213 (Fax)

Recommended Sterilization Procedures for Zeus® Intervertebral Body Fusion Instrumentation and Implants Provided Non-Sterile:

**Manufacturer:** Amendia, Inc.

**Method:** Manual Cleaning and Steam Sterilization

**Device(s):** Trays/Implants/Instruments

Cautions:	<p>The Zeus® Intervertebral Body Fusion systems provided NON-STERILE should be cleaned and sterilized before use.</p> <p>Automated cleaning may not be effective. A thorough, manual cleaning process is recommended.</p> <p>Cleaning agents with chlorine or chloride as an active ingredient are corrosive to stainless steel and should not be used.</p> <p>Saline solution has a corrosive effect on stainless steel and should not be used.</p> <p>Use only neutral pH cleaning agents and detergents.</p> <p>All Zeus® Intervertebral Body Fusion IMPLANTS are single use. Therefore these guidelines are not intended for USED Zeus® spinal implants or DISPOSABLE, single use instruments.</p> <p>The Zeus® Intervertebral Body Fusion System has not been evaluated for safety and compatibility in the MR environment. The Zeus® Intervertebral Body Fusion System has not been tested for heating or migration in the MR environment.</p>
Limitations on Reprocessing:	<p>Repeated processing has limited effect on REUSABLE instruments.</p> <p>End of life is normally determined by wear and damage due to use.</p>
Instructions:	
Point of Use:	Use clean flowing water and disposable wipes to remove excess soil. Reprocess instruments as soon as possible to prevent body fluid and tissue from drying on instruments prior to cleaning.
Preparation for decontamination:	Disassemble all components to provide maximum exposure for cleaning.
Cleaning - Automated:	Automated washer/disinfector systems are not recommended as the sole cleaning method for surgical instruments. An automated system may be used as a follow-up method to manual cleaning.
Cleaning- Manual:	<ol style="list-style-type: none"> <li>1. Disassemble all components before cleaning.</li> <li>2. Completely submerge instruments in enzyme solution and allow to soak for a minimum of 20 minutes. Use a soft-bristled, nylon brush to gently scrub the device until all visible soil has been removed. Particular attention must be given to crevices, lumens, mated surfaces, connectors and other hard-to-clean areas. Lumens should be cleaned with a long, narrow, and appropriately sized soft-bristled brush (e.g. pipe cleaner brush).</li> <li>3. Remove the devices from the enzyme solution and rinse in tap water for a minimum of 3 minutes. Thoroughly and aggressively flush lumens, holes and other difficult-to-reach areas.</li> <li>4. After manual cleaning, and all visible blood, soft tissue, and bone have been removed ultra-sonic cleaning may be used. Place prepared cleaning agents in a sonication unit. Completely submerge device in cleaning solution and sonicate for a minimum of 10 minutes at 45-50kHz.</li> <li>5. Rinse instrument in purified water for at least 3 minutes or until there is no sign of blood or soil on the device or in the rinse stream. Thoroughly and aggressively flush lumens, holes and other difficult-to-reach areas. Use de-ionized water for final rinse of all components.</li> </ol>

# Important Product Information

Cleaning- Manual Continued...	<p>6. Repeat the sonication and rinse steps above until all visible contamination has been removed.</p> <p>7. Thoroughly and promptly, remove excess moisture from the instrument with a clean, absorbent and non-shedding wipe. Allow the tray and components to dry for a minimum of 15 minutes. The tray and components must be thoroughly dry prior to sterilization cycle.</p>
Disinfection:	Disinfection is only acceptable as an adjunct to full sterilization for reusable surgical instruments.
Maintenance, inspection, and testing;	<p>Carefully inspect each device to ensure that all visible blood and soil have been removed.</p> <p>Inspect lumens to confirm that all foreign material has been removed.</p> <p>Visually inspect for damage and/or wear.</p> <p>Note: If any damage or wear is noted that impairs the function of the instrument, contact your Amendia representative for a replacement.</p>
Packaging:	This set of components may be loaded into a dedicated tray, supplied by the manufacturer, for sterilization.
Sterilization:	<p>Visually inspect all components for any remaining debris prior to sterilization.</p> <p>The Zeus® Intervertebral Body Fusion system components provided NON-STERILE should be autoclave sterilized using the sterilizer manufacturer's instructions and the institution's procedures for ensuring sterility. The sterilization cycle should occur in a calibrated autoclave.</p> <p>Zeus-A® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>Zeus-C® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>Zeus-T® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>Zeus-P® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>Zeus-V® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>Zeus-L® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 10 minutes at 270°F (132°C.)</p> <p>Zeus-O® system components should be sterilized utilizing a pre-vacuum steam autoclave for a minimum of 4 minutes at 270°F (132°C.)</p> <p>The 10 minute, 270° pre-vacuum steam sterilization cycle is not considered by the Food and Drug Administration to be a standard sterilization cycle. It is the end user's responsibility to use only sterilizers and accessories (such as sterilization wraps, sterilization pouches, chemical indicators, biological indicators, and sterilization cassettes) that have been cleared by the Food and Drug Administration for the selected sterilization cycle specifications (time and temperature).</p>
Drying:	<p>A minimum drying time of 20 minutes, after sterilization, is recommended.</p> <p>Drying times may vary according to load size and should be increased for large loads.</p> <p>Dry, thoroughly and promptly, after both cleaning and sterilization.</p>
Storage:	Store components in a clean, dry, limited access area that is well ventilated and provides protection from dust, moisture, insects, vermin, and extremes in humidity and temperature.

The instructions provided above have been validated by Amendia as being CAPABLE of preparing a medical device for reuse. It remains the responsibility of the processor to ensure that the reprocessing as actually performed, using equipment, materials, and personnel in the reprocessing facility, achieves the desired result. This normally requires validation and routine monitoring of the process. Any deviation by the re-processor from these instructions should be properly evaluated for effectiveness to avoid potential adverse consequences.



1755 West Oak Parkway  
Marietta, GA 30062  
Phone: 877-755-3329  
Fax: 877-420-1213  
[www.Amendia.com](http://www.Amendia.com)  
[info@amendia.com](mailto:info@amendia.com)

510(k) Number: K081614  
MM-025, Rev 1