




SELF-ADHERENT SOFT SILICONE DRESSING FOR SCAR CARE

- No extra fixation is needed
- Conforms well to body contours
- Can be worn in daily activities

Film/nonwoven backing

- Moisture vapour permeable
- Waterproof

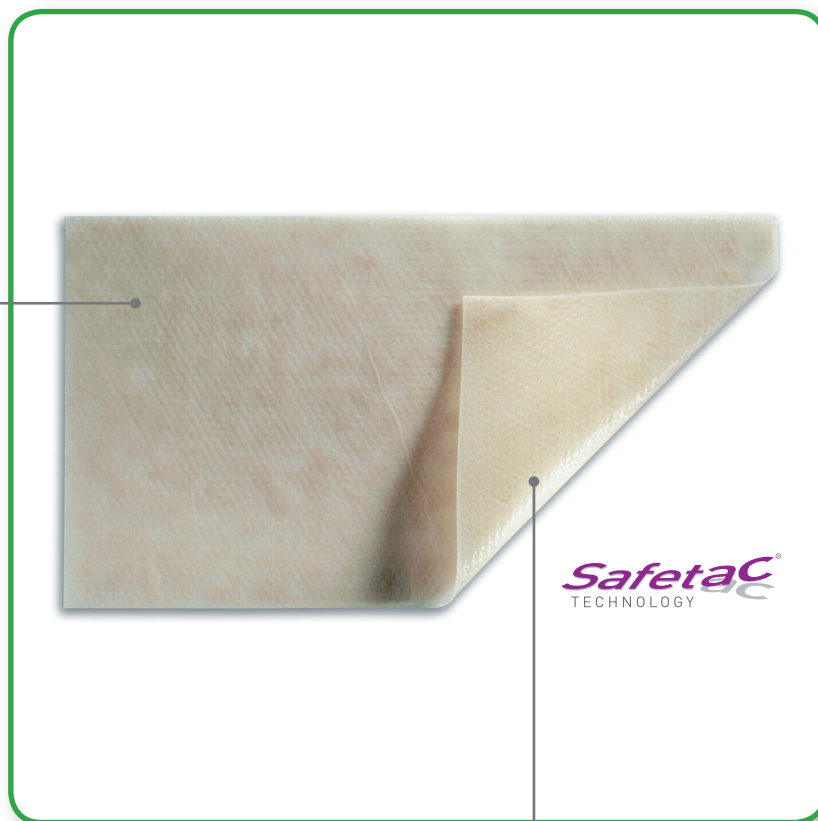


Safetac® layer (soft silicone)

- Atraumatic to the skin on removal


HOW MEPIFORM WORKS

Topical silicone has empirically been shown to have positive impact on hypertrophic scars and keloids.¹ It may take from 3 months up to a year or more to improve an old scar, depending on the condition of the scar tissue. Mepiform can optimally be worn 24 hours a day. It is recommended that Mepiform is removed once a day, e.g. when showering or bathing, for inspection and washing of the skin. For prophylactic treatment, Mepiform can be used for 2–6 months depending on the condition of the scar.


HOW TO USE MEPIFORM

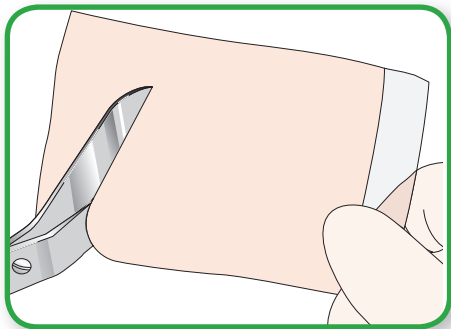

1. Remove the release film. If needed, cut to shape.release film.

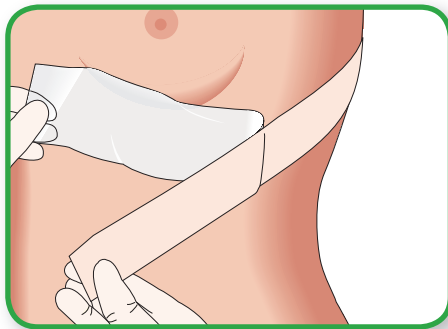


2. Apply Mepiform to the scar.

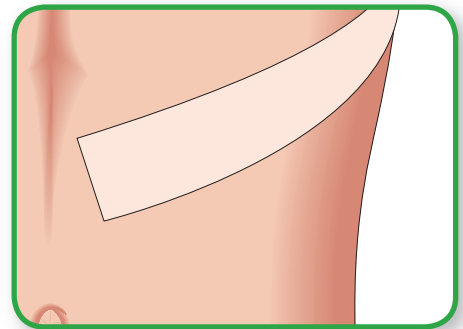

3. Reapply after daily cleaning.

BENEFITS OF MEPIFORM

- No extra fixation is needed
- Mepiform is thin, flexible and discreet
- Can be worn in daily activities
- Can be reapplied
- One Mepiform product can be used for several days

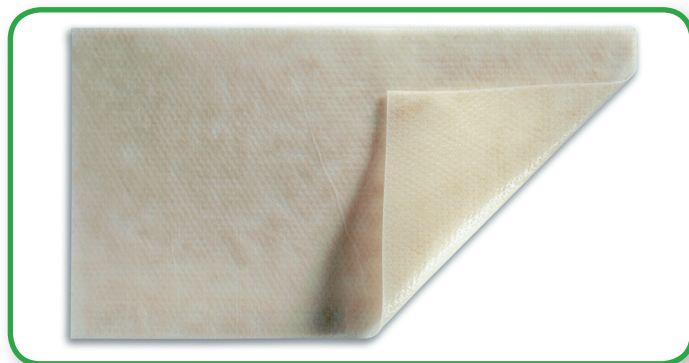

NOTE

Should maceration or rash occur, allow the skin to rest until the symptom has disappeared. Continue treatment and gradually increase the treatment time per day. Mepiform has an ultraviolet protection factor equivalent to 7.7.²

MEPIFORM ORDERING INFORMATION

Art. No	Size cm	Pieces Per Inner	PIP Code	NHS Code
293250	5 x 7	5	332-7335	ELY200
293450	9 x 18	5	332-7343	ELY201
293150	4 x 31	5	332-7350	ELY199

References:

1. Quinn KJ (1987), Burns 13:33-40; Ahn ST et al (1991), Arch Surg 126:499-504; Dockery GL et al [1994, J Foot Ankle Surg 33:110-19; Katz BE(1995) Cutis 56:65-7; Carney SA et al [1994] Burns 20:163-7.

2. AS/NZS 4399:1996.



Mölnlycke Health Care Ltd, Arenson Centre, Arenson Way, Dunstable, Bedfordshire LU5 5UL, UK. Tel. +44 (0) 870 6060766, www.molnlycke.co.uk

The symbol and the word mark are both registered trademarks or trademark pendings of Mölnlycke Health Care. Clinical references available on request.

UK175560611