

Safe Patient Handling Solutions

McAuley Medical, Inc.

THE **SOFT GOLD ROLLBOARD** Patent Pending

HEAT-SEALED COVER – NO NEEDLE AND THREAD STITCHING

*Don't be misled by imitators -
insist on these features:*

Heat-Sealed Cover • Anti-Microbial
3-Year Warranty • Anti-Static • MRI Safe
Molded Plastic Handles



1.



2.

See short in-service videos on our website.

Available with
**FREE In-Service
Video and
e-Learning
Module**



3.

Rollboard bridges the "Gap"

Older airflow technology requires a bridge to cover gaps between surfaces. Slide sheet transfers across gaps are painful to the patient. Rollboard automatically bridges the gap and the patient feels nothing.

Environmentally Friendly

- NO Electricity • NO Noise
- NO Laundry/Detergents required



**ALL COVERS, SMOOTH PLASTIC HANDLES
AND INNER CORE USE STATE-OF-THE-ART
ANTI-MICROBIAL TECHNOLOGIES**

Gold Fabric:

- *Ultra-fresh*® additive helps keep our gold fabrics fresh, hygienic and odor-free.



Plastic Handles:

- *Biomaster*® anti-microbial additive provides a slow release of silver ions, safely inhibiting bacterial growth.



Inner Core:

- Molded *ARPRO*® meets the following ASTM fungus, mildew and bacteria resistance specifications: C1338, G-21 and G-22.



Our unique "Conveyor Belt" action...



carries the patient from
surface to surface,
powered only by a
gentle "push"



Made with pride
in the USA

Soft Foam Core - NO rollers

Lateral transfer *without lifting* No lifting required, staff back injury is reduced. The result can be savings on both worker's compensation and insurance.

Lateral transfer *without noise* Patients are not disturbed, annoyed or frightened by noise associated with airflow systems.

Priced well below other technologies You can purchase a Gold Rollboard for every room in your facility at a fraction of the cost of mechanical devices and airflow mattresses.

Patient rides on a soft surface Patient is comfortable with far less risk of skin tears, injury, bruising and agitation associated with transfers on slide sheets and plastic boards.



NiN Healthcare Kaiser Permanente Contract # KP81901
Please place your order via one link
OR email orders@ninhealthcare.com

McAuley Medical, Inc.

888.399.7287 • 352.399.6816 • FAX: 603.227.7286 • WWW.MCAULEYMEDICAL.COM

The complete family of exclusive McAuley Medical **SOFT Gold Rollboards**

Gold Short OR
30" x 15"
Part # 5430



Gold OR
48" x 15"
Part # 5400



Gold Birthing
34" x 19"
Part # 5300



Gold Bari
68" x 19" Part # 5200



Gold Standard
68" x 19"
Part # 5100



Gold Radiology
68" x 15"
Part # 5500



Gold EMS
50" x 15"
Part # 5900





Gold Rollboard Disposable Cover Sheet

For medium duty transfers, McAuley Medical offers this lower-price-point Blue Board.

Short Wide
Part # 8304
34" x 19"



Long Wide
Part # 8300
68" x 19"



Short Narrow
Part # 8301
48" x 15"



Long Narrow
Part # 8305
68" x 15"





For lighter duty applications, we offer the Banana ROLASLIDE.

Part # 4420
69.5" x 21.5"



The McAuley Rollboard Collection features:


	Blue Rollboard	Gold Rollboard
Heat-Sealed Cover	✓	✓
Smooth, Anti-Microbial Handles	✓	✓
Anti-Static cover	✓	✓
Can also be used in Fowler positions	✓	✓
Molded plastic handles	✓	✓
Heavy duty use	✓	✓
MRI Compatible	✓	✓
Anti-Microbial Core	✓	✓
Anti-Microbial cover		✓
Models Available	4	7
Warranty	1-year	3-year



Bendable Patient Transfer Board
Anti-static, low-friction plastic allows for easier lateral transfer.
Part # 7001

GOLD ROLLBOARD MODELS ARE DESIGNED TO FIT A VARIETY OF APPLICATIONS THROUGHOUT YOUR FACILITY.

Patent Pending

Gold Rollboard Model	Application(s) (Recommendations)	Department(s)
A. GOLD OR 48" x 15" Part # 5400	To and from OR Table (1 for each OR)	Operating Room
B. GOLD OR SHORT 30" x 15" Part # 5430	To and from OR Table (1 for each OR)	Operating Room
C. GOLD Birthing 34" x 19" Part # 5300	To and from Bed, Stretcher and OR Table (1 for each 6 birthing area beds, if short board is preferred)	Birthing Unit
D. GOLD Bari 68" x 19" Part # 5200	To and from Bariatric Bed and Stretcher (reinforced to carry more weight) (1 for each Bariatric Suite)	ED, ICU, Imaging, Inpatient Floor, LTC, Homecare
 The following boards are foldable to facilitate Fowler Position Transfers:		
E. GOLD Standard 68" x 19" Part # 5100	To and from Bed, Stretcher and XRay Table (1 for each 8 inpatient beds and 8 ICU beds) (1 for each XRay room and morgue)	ED, ICU, Imaging, Birthing, Inpatient Floor, LTC, Homecare, Ultrasound, Morgue
F. GOLD Radiology 68" x 15" Part # 5500	To and from MRI table, CT, Cath Lab, OR Table (1 for each MRI, CT and Cath Lab table)	MRI, CT, Cath Lab, IR, Nuclear Medicine, OR
G. GOLD EMS 50" x 15" Part # 5900	EMS Stretcher to Bed (1 per 6 ER beds and 1 per Ambulance)	Ambulance, ED

LATERAL TRANSFER - ROLLBOARDS



MRI SAFE!
No air holes,
No Velcro

No place for microbes to hide.

Hangs on wall hooks

System is easy to store, takes up little space, and can hang conveniently in every room of your facility.

Also available:

Disposable Cover Sheets

Long Short
Part # 1000 Part # 1400
83" x 23" 63" x 23"



Anti-Static Material

Our Gold EMS Rollboard was designed specifically to meet the needs of Emergency Response Professionals, reducing back injuries AND workers comp costs.



New Handles: All McAuley Rollboards feature sleek, smooth plastic handles that are easy to grip, easy to clean, and release silver anti-microbial ions. They represent a vast improvement over the handle designs shown below.



Handy Rollboard Accessories

Disposable cover sheet



Inexpensive, disposable thin plastic cover protects Rollboard from soiling and provides an additional layer of infection control.

Comes in rolls of 50; two sizes available:

Short, part # 1400, 63" x 23"

Long, part # 1000, 83" x 23"




NiN Healthcare Kaiser Permanente Contract # KP81901
Please place your order via one link
OR email orders@ninhealthcare.com

Disposable Cover Dispenser and Hanger System




Rollboard Disposable Cover Dispenser
Part # 1088

- FDA Approved Powder Coating • Steel construction
- Can sit on a counter top (rubber feet included)
- Easy to wipe clean

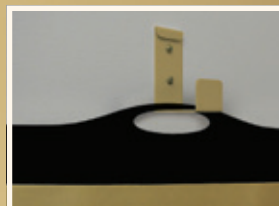

3" Long Wall Mount Hook
Part # 1086

MRI Safe Wall Mount Hook
Part # 2086

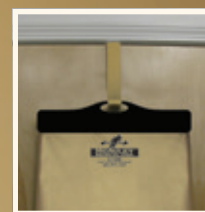

8" Long Door Mount Hanger
Part # 1087

MRI Safe Door Mount Hanger
Part # 2087

20" Long Door Mount Hanger
Part # 1089

Closet Rod Hanger
Part # 9999

THE SCIENTIFIC EVIDENCE - TWO MAJOR STUDIES:

RIGHT: MEMIC, a leading New England Worker's Compensation insurance specialist, recently completed a rigorous scientific evaluation of Lateral Transfer Devices to find the best solutions for reducing friction and lowering push/pull forces, and a review of non-motorized friction reducing devices used for Lateral Transfer. The results speak for themselves:

BELOW: Here are the results of a recent study by Andrea Baptiste, MA, OT, CIE, and Kay Steadman, MA, OTR, CHSP.

The surprising fact is that in three out of four of the categories, the Rollboard not only required the least force, but also required as much as 15% less force than the Air Assisted product!

Once again the Rollboard proves to be the safest Lateral Transfer Device and, at about a quarter of the cost of air-assisted devices, the best value for your safe patient handling dollar!

For the complete details on these two important studies, visit www.mcauleymedical.com

Lateral Transfer Devices that are not motorized

Non Motorized Device	Peak Pull Force Average of 4 trials	Patient Comfort Rating 1=best	Probability of Low Back Injury
Cotton Draw Sheet	Highest Peak Pull Force		Highest Probability of Back Injury
Wooden Board			
White Plastic Board		4	
Blue Plastic Transglide Board		2	
Metal Roller Board		7	
Plastic Slide Sheet		6	
Two Piece Nylon Slide Sheet		5	
Tubular Nylon Slide Sheet	Lowest Peak Pull Force	3	Lowest Probability of Back Injury
Rollboard		1	

All of these devices cost less than \$1,000 each

Ranked in order of *least* (1) to *most* (6) physically stressful or demanding.

Peak force when transferring 171 lbs

	Device	TASK Duration (seconds)	Peak Force (Resultant) Newtons	Rank
S H O U L D E R	Incontinent Pad/draw sheet	9.21	259.22	6
	Incontinent pad over slide board	5.63	189.04	5
	Friction reducing sheet	5.38	183.58	3
	Slide board	6.03	180.03	4
	Rollboard	5.43	154.10	1
	Air assisted	5.28	177.20	2

Peak force when transferring 223 lbs

	Device	TASK Duration (seconds)	Peak Force (Resultant) Newtons	Rank
S H O U L D E R	Incontinent Pad/draw sheet	9.21	284.98	6
	Incontinent pad over slide board	5.63	208.87	3
	Friction reducing sheet	5.38	247.25	5
	Slide board	6.03	208.68	2
	Rollboard	5.43	199.09	1
	Air assisted	5.28	218.90	4

Peak force when transferring 223 lbs

	Device	TASK Duration (seconds)	Peak Force (Resultant) Newtons	Rank
L O W B A C K	Incontinent Pad/draw sheet	9.21	724.76	6
	Incontinent pad over slide board	5.63	623.83	4
	Friction reducing sheet	5.38	674.87	5
	Slide board	6.03	618.85	3
	Rollboard	5.43	604.63	1
	Air assisted	5.28	615.12	2