



UNOLOCK® INFUSION PRODUCTS

Now, Safety and Performance are Easy

UNOLOCK® Infusion Sets include a line of both conventional and safety infusion sets and feature siliconized needles, 100% vision camera inspection, and color coded wings. UNOLOCK® Sets are manufactured to strict quality standards and contain unique features to improve user convenience, safety, and performance.



Bringing **Reliability, Innovation and Value** to Healthcare®

UNOLOK® Infusion Sets


The Performance and Convenience You Need

- Siliconized short bevel needle provides smoother insertion
- Computerized grinding and 100% vision camera inspection ensures sharp, consistent, high quality needles
- Flexible, color coded wing for improved tactile feel and easy gauge identification
- Safety Sets include single handed safety shield activation for convenience and ease of use
- Manufactured in FDA inspected facilities and conform to ISO 9262, 594, and 13485

A Range of Sizes To Meet Your Needs

- Sizes: 19, 21, 23, and 25 gauges. Safety Sets are 8", Conventional Sets are 12" (Sterile)

UNOLOK® Infusion Sets Gauge Color Chart

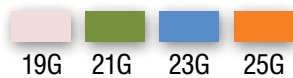

UNOLOK® Safety Infusion Set:



1. Prior to needle removal from a patient, the safety shield should be positioned so it is touching the back of the wings
2. As the needle is removed from the patient, the shield should be slid forward over the wings. When it **Clicks**, its safe

UNOLOK® Conventional Infusion Set:



To Inquire about MYCO's Products or receive complimentary evaluation samples, please contact:

 **1.800.454.6926 or 1.919.460.2535**
 **Info@mycomedical.com**



Bringing **Reliability, Innovation and Value** to Healthcare®