


# VIO® 100 C



**A COMPACT ERBE  
GENERATOR for  
less-invasive procedures.**

## FEATURES:

- ✦ Compact/Lightweight
- ✦ Up to 4 storable programs
- ✦ Return Electrode Monitoring (REM)

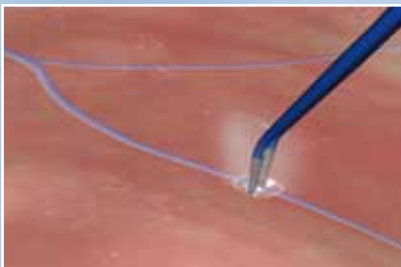

## MODES AVAILABLE:

- ✦ AUTO CUT® – Electrosurgery (ES) cutting with minimal to medium hemostasis.
- ✦ DRY CUT – ES cutting with intense hemostasis.
- ✦ SOFT COAG – Non-sparking, soft monopolar contact coagulation with minimal sticking and carbonization.
- ✦ FORCED COAG – Effective, fast, coagulation.
- ✦ BIPOLAR SOFT COAG – Non-sparking, soft bipolar contact coagulation with minimal sticking and carbonization.

## PROPERTIES:

- ✦ BIPOLAR COAGULATION
- ✦ MONOPOLAR CUTTING
- ✦ MONOPOLAR COAGULATION

**ERBE**  
USA INCORPORATED  
Surgical Systems

# ERBE VIO® 100 C

## Product Data

### DESCRIPTION

### PART NUMBER

ERBE VIO® 100 C Electrosurgical Unit, 120 V / 60 Hz, UL-Version

10140-500



VIO C Two (Double/Dual) Pedal Footswitch, AP and IP X8 Equipment

20189-107



\* Also Available: VIO C One (Single) COAG Pedal Footswitch, AP and IP X8 Equipment

20188-102

## Technical Data

### POWER OUTPUT

**Maximum CUT output** Up to 100 Watts

**Maximum COAG output** Up to 80 Watts

### MAIN CONNECTION

**System Voltage** 100-120 / 220-240 V

**Power Frequency** 50 / 60 Hz

### DIMENSIONS AND WEIGHT

**Width x Height x Depth** 11" x 5.3" x 11.8"

**Weight** 8 lbs 13 oz

### STANDARDS

**Classification according to EC Directive 93/42/EEC** IIb

**Protective class according to EN 60 601-1** I

**Type according to EN 60 601-1** CF

**ERBE**  
USA INCORPORATED  
Surgical Systems

2225 Northwest Parkway  
Marietta, GA 30067

Tel. 800.778.3723  
Tel. 770.955.4400  
Fax. 770.955.2577  
[www.erbe-usa.com](http://www.erbe-usa.com)

MKT/5057/00 (02/11)