


Eclipse[®]

Super Absorbent Dressing Range

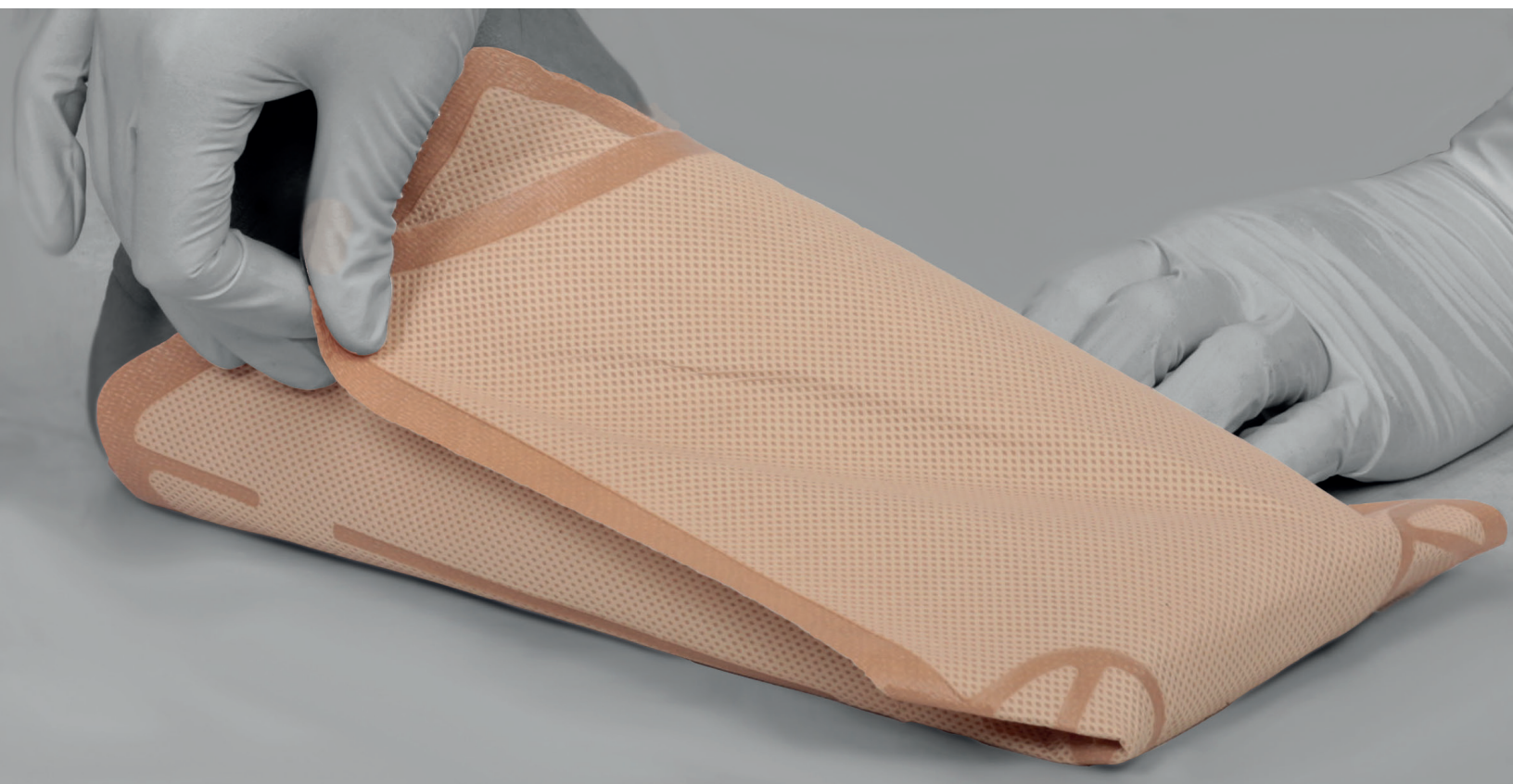

Effective exudate management

Distributed by




PASSION. PARTNERSHIP. POSSIBILITIES.

www.dukal.com



www.advancismedical.com

How do the layers of the Eclipse® dressings work?

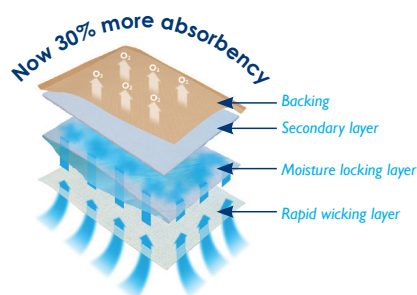

The backing

The water-resistant barrier backing on Eclipse® dressings are bacteria and viral proof and critically prevent strike-through.



Secondary layer

This layer sits underneath the backing and on top of the moisture locking layer.



Moisture locking layer


The central layer of the Eclipse® dressing consists of a sheet of highly absorbent crystals and a mechanically bonded cellulose pad.



Rapid wicking layer

This layer draws up the exudate from the wound bed delivering the exudate into the moisture lock layer.

Manufactured in our state-of-the-art facility in the UK, our range of Eclipse® super absorbent dressings deliver quality performance in the controlled management of exudate. Our recent developments in the Eclipse® range of dressings have seen an increase in the absorbency of our products by 30%, with no impact on cost.

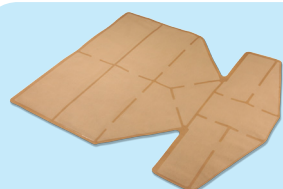

Eclipse Boot®

Super absorbent leg wrap dressing

A unique clinical solution to assist in the management of excess fluid often associated with lower leg oedema.

Eclipse Boot® is a super absorbent dressing specifically designed to manage high levels of fluid in the lower limb.

The preshaped design overcomes the traditional time-wasting complexities of dressing 'leaky legs' with a combination of dressings that are often bulky and restrict patient mobility.

- ✓ Anatomically designed to fit the lower limb and foot
- ✓ 12 Super absorbent interconnected compartments to optimize absorbency
- ✓ Highly conformable
- ✓ Significantly improves patient mobility
- ✓ Easy to apply, single dressing design
- ✓ Reduce dressing application time
- ✓ Available in three sizes



Eclipse Foot®

Super absorbent foot wrap dressing

A super absorbent dressing, designed to be applied securely around the foot area to manage wound exudate.

The dressing utilizes the high fluid handling capabilities of the Eclipse® range in a preshaped dressing designed to fit comfortably and effectively around the foot. It is simple to apply and offers enhanced patient comfort.

- ✓ Rapid absorption
- ✓ High capacity
- ✓ Fluid repellent backing
- ✓ Strength and durability
- ✓ Stay dry technology
- ✓ Thin and conforming
- ✓ Particularly useful in issues related to the treatment of the diabetic foot



Item #	Description	Packaging	HCPCS
CR4014	Eclipse Boot®, 24" x 28"	5/bx	A6198
CR4219	Eclipse Foot®, 13" x 19"	5/bx	A6198
CR4229	Eclipse Boot®, 28" x 32"	5/bx	A6198
CR4237	Eclipse Boot®, 22" x 19"	5/bx	A6198

Your patient outcomes can be improved when Eclipse® is used in conjunction with our range of 100% Medical Grade Manuka Honey products.



DUKAL Corporation
2 Fleetwood Court
Ronkonkoma, NY 11779
Tel: 1-800-243-0741 | Email: advancisinfo@dukal.com

MAR434

Advancis
Medical
www.advancismedical.com