


# Adjustable height staples. **Wide range of tissue.**

## PROXIMATE® Reloadable Linear Stapler (TL)

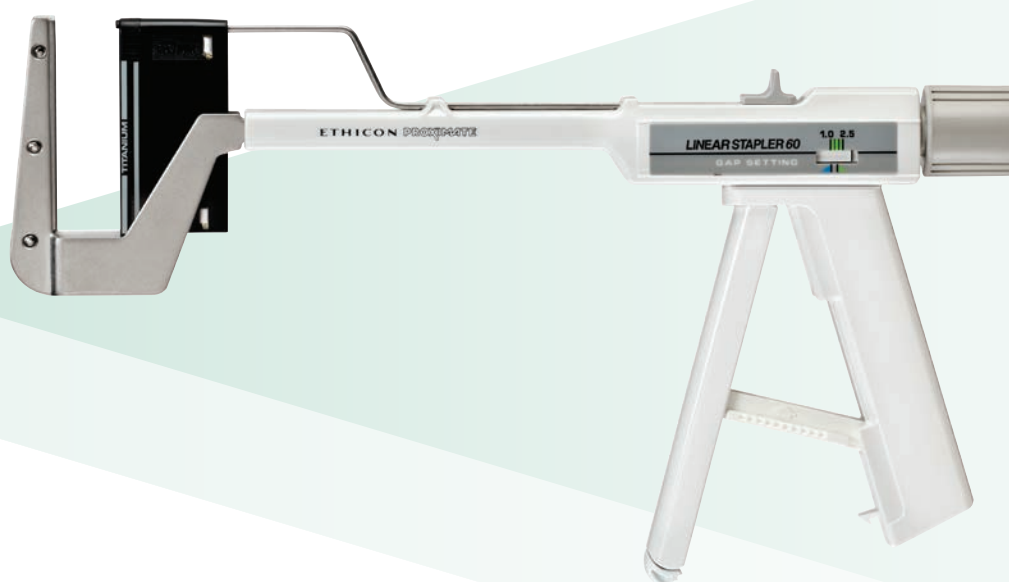

### Consistent staple formation

- Parallel jaw closure helps to compress tissue evenly within the jaws of the instrument for consistent staple formation

### Adjustable staple height

- Adjusting Knob allows surgeon to choose the staple height to better accommodate various tissue thicknesses

### Applicable specialties

- General
- Colorectal
- Thoracic

## Product Overview

- Heavy wire option provides the ability to staple thick tissue
- Two staggered rows of titanium alloy staples in 30mm, 60mm, 90mm sizes

## Device Specifications

Description	Device code	Tissue	Staple rows	Staple line length	Number of staples	Max firings	Qty/box
Proximate Linear Stapler 30mm	TL30 TLH30	Regular/Very thick	2	30mm	11	4	3
Proximate Linear Stapler 60mm	TL60 TLH60	Regular/Very thick	2	60mm	21	4	3
Proximate Linear Stapler 90mm	TL90 TLH90	Regular/Very thick	2	90mm	33	4	3



## Reload Specifications

Description	Reload code	Tissue	Closed staple height	Wire diameter	Open staple height	Qty/box
Proximate Linear Stapler 30mm Reload	TR30 TRH30	Regular/Very thick	1.0-2.5mm	0.23mm 0.30mm	4.4mm 4.8mm	12
Proximate Linear Stapler 60mm Reload	TR60 TRH60	Regular/Very thick	1.0-2.5mm 1.5-2.5mm	0.23mm 0.27mm	4.4mm 5.5mm	12
Proximate Linear Stapler 90mm Reload	TR90 TRH90	Regular/Very thick	1.0-2.5mm 1.5-2.5mm	0.23mm 0.27mm	4.4mm 5.5mm	12

