


### The Stonehaven Medical BAL2000 Pacific 6 Section Hi-Lo Treatment Table

is designed with 3 actuators or motors for complete control of cushion articulation, as well as table height adjustment. Ease of operation with Under Cushion Finger Touch Controls and 4 Elevating Bars for Table height adjustment, the Stonehaven BAL2000 Pacific 6 Section is the most effortless free operative Treatment Table available. The Head and Back section adjust from -28 to 80 degrees and the center section from 0 to 20 degrees. The Arm Rest, 2, adjust outward as well as up and down for patient comfort. Constructed of heavy-duty steel frame with lifting capacity of 500 lbs. Cushions are 2" foam and choice of 6 attractive upholstery colors. The Head Section includes Breathing hole and plug. Dual Foot operative Casters, 4, are a standard feature.

### SPECIFICATIONS

- Gas Struts, hand lever activated, providing smooth adjustment at movable top sections
- Cushion width of 28" (71 cm) x length of 81" (206 cm)
- Dual foot lever retractable heavy-duty hospital grade casters (4)
- FDA/CE listed with hospital grade electrical cord
- ADA Compliant (see IRS, Section 44 or tax consultant regarding tax credit allowance)
- Frame designed for patient assisted lift accessibility, for ease of patient transfer
- Cushions secured to the heavy-duty metal frame with an integral T-Nut Bolt assembly
- Hands free, foot operated up/down elevation levers
- Complete with Anti Slipping Foot levelers

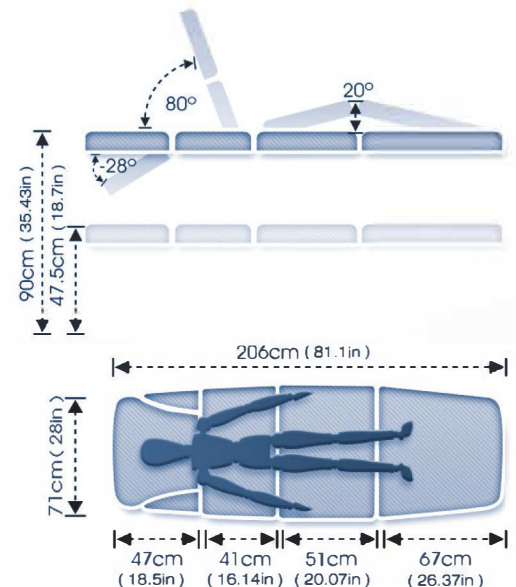

### OPTIONS

- Custom Upholstery width and color (CU)
- Foam Memory Face Pad (FMP)
- Hydraulic Elevation Pump (HP)
- Paper Dispenser (PD)
- Foot or Hand Control Switch, in lieu of Elevating Bars



\*note special order

We reserve the right to modify designs without prior notice

**Stonehaven Medical, Inc.**

22607 La Palma Ave., Suite 408, Yorba Linda, CA 92887

P: 714-692-9200 • 888-215-4205 • F: 714-692-9800

[www.stonehavenmedical.com](http://www.stonehavenmedical.com)