

AT4


Tourniquet System

Ideal for single or bi-lateral
orthopaedic surgery and pain
management



AT4

Pneumatic
Tourniquet

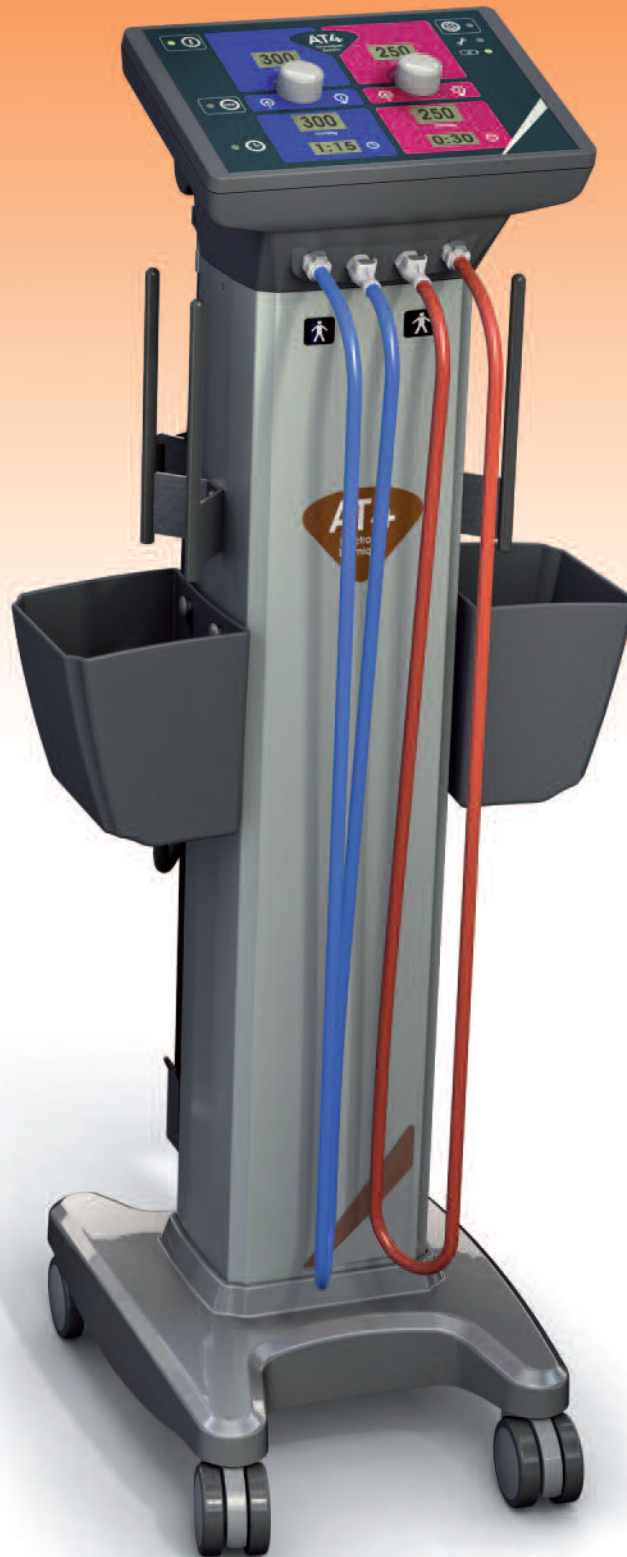

The AT4 Pneumatic Tourniquet:

A traditional dual channel pneumatic tourniquet operating from 4 or 7 bar piped air service comes with mobile stand and cuff supply hoses

NEW

NEW

AT4⁺

Electronic
Tourniquet



The AT4 Electronic Tourniquet:
A state-of-the-art dual channel electronic
touniquet with integral power supply which requires
no external compressed air supply comes with
mobile stand and cuff supply hoses

Innovative technology – practically applied

The AT4 Tourniquet System is a completely new design which is intuitive and simple to use. Its compact construction is light weight yet robust, containing advanced technology which offers maximum safety and reliability.

The electronic version operates from its own integral air supply, removing the need for an external pressure source, and both this and the pneumatic version have sensitive leak compensation facilities as well as automatic self-diagnostic checks on start-up, which means they are calibrated every time they are switched on.

An innovative added safety feature has been incorporated for Bier's block regional anaesthesia procedures (IVRA): it will simply not be possible to deflate cuffs in the wrong order risking the dangerous release of anesthetic into the rest of the patient's body.



Automatic self diagnostic checks
On start-up and throughout use



Independent displays
Pre-select and actual cuff
pressure readings



Precision controlled regulators
User friendly simple to use controls



IVRA (Bier's block)
Procedural guidance feature



Automatic LCD timers
With user defined elapsed time alarm



Additional storage facility



Cuff hooks
For ease of storage



For more information on any of our products or service contracts, call:

1800.499.5370

www.dremed.com



QUALITY ASSURED



GLOBAL DISTRIBUTION

The AT4 Electronic Tourniquet:


A state-of-the-art dual channel electronic tourniquet with integral power supply which requires no external compressed air supply comes with mobile stand and cuff supply hoses

The AT4 Pneumatic Tourniquet:

A traditional dual channel pneumatic tourniquet operating from 4 or 7 bar piped air service comes with mobile stand and cuff supply hoses



See back cover for optional extras

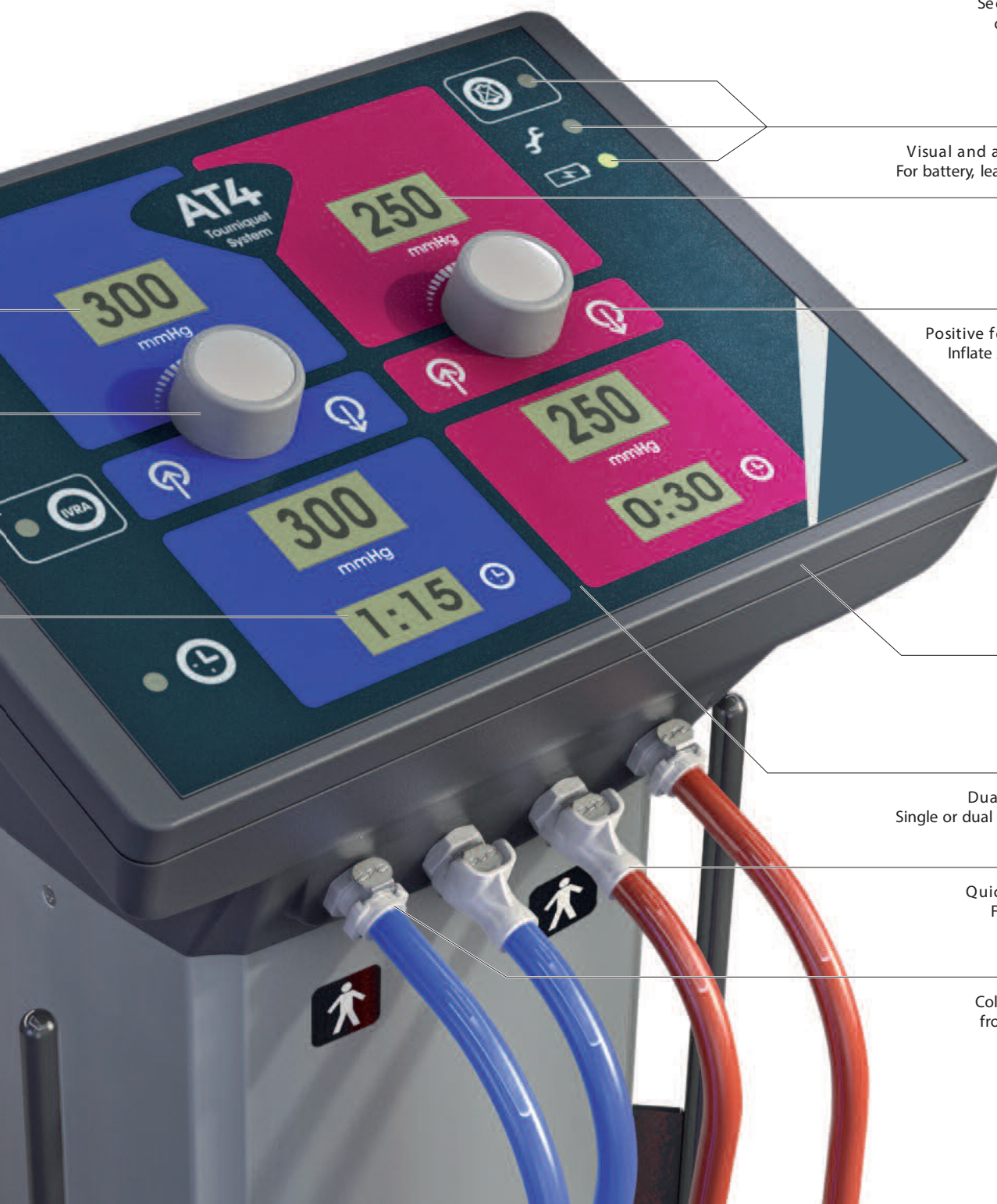

Visual and audible warnings
For battery, leaks and calibration



Pressure range
0 – 600 mmHg



Positive feel push buttons
Inflate / deflate functions



Easy to clean
Low maintenance



Dual channel system
Single or dual channel operation



Quick release hosing
For ease of storage



Color coded channels
from controls to cuff

