

RAPID DIAGNOSTICS

[www.osomtests.com](http://www.osomtests.com)

# osom<sup>®</sup>

## Trichomonas Rapid Test

FOR THE DETECTION OF *TRICHOMONAS VAGINALIS* ANTIGENS



### FEATURES:

- ✓ **ACCURATE;** 83% sensitivity vs. culture; 95% agreement with culture and wet mount combined
- ✓ **FAST;** Results in 10 minutes or less
- ✓ **SAMPLE;** Vaginal Swabs
- ✓ **OBJECTIVE;** Detects the antigen; removes uncertainty associated with live organism instability
- ✓ **TRAINING;** Online training and certification; [www.osomtraining.com](http://www.osomtraining.com)
- ✓ **CUSTOM CARE;** Technical support staffed with certified medical technologists

## RAPID DIAGNOSTICS

[www.osomtests.com](http://www.osomtests.com)


# osom<sup>®</sup>

## Trichomonas Rapid Test

**Sekisui Diagnostics, LLC (AMERICAS)**  
4 Hartwell Place Lexington, MA 02421  
Phone: 800 332 1042  
Fax: 800 762 6311  
Email: [osomrapidtests@sekisui-dx.com](mailto:osomrapidtests@sekisui-dx.com)  
ISO 13485 Certified

**Sekisui Diagnostics (UK) Limited (INTERNATIONAL)**  
Liphook Way, Allington, Maidstone, Kent ME16 0LQ, UK  
Phone: +44 1622 607800 Fax: +44 1622 607801  
Email: [info@sekisui-dx.com](mailto:info@sekisui-dx.com)  
ISO 9001 Certified, ISO 13485 Certified

### TEST PROCEDURE / RESULTS

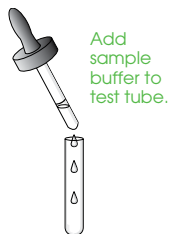

Insert swab and mix.



Let stand 1 MINUTE.



Insert test stick into extracted sample.



Read results at 10 MINUTES.



#### DESCRIPTION

OSOM<sup>®</sup> TRICHOMONAS RAPID TEST  
TRICHOMONAS POSITIVE CONTROL KIT

OSOM<sup>®</sup> VAGINITIS TRANSPORT SYSTEM (TRICH & BV)

#### CONFIGURATION

25 TESTS

50 SWAB SETS

#### CATALOG #

181 OR 181E (CE MARKED)  
182 (CE MARKED)

185

#### STORAGE

ROOM TEMP  
ROOM TEMP

ROOM TEMP

#### CPT CODE

87808QW

