

800-525-9339 • 651-674-8301
www.swedeo.com

- Easily control support with external figure eight straps
- Stays tight and secure with exclusive lacing system and top strap
- Cool, comfortable, breathable tongue
- Durable, lightweight ballistic nylon is ideal for sports
- Full elastic back reduces irritation and potential for blistering
- Bilateral design fits both right and left ankle
- Industry leading 12 month warranty, satisfaction guarantee



Sized by US shoe size

BLACK MODEL #	WHITE MODEL #	SIZE	MEN'S	WOMEN'S
AKL-6334-BK-1XS	AKL-6334-WH-1XS	XS	3-5 1/2	4-6 1/2
AKL-6334-BK-SML	AKL-6334-WH-SML	SM	6 - 7 1/2	7-9 1/2
AKL-6334-BK-MED	AKL-6334-WH-MED	MED	8-10 1/2	10-11 1/2
AKL-6334-BK-LRG	AKL-6334-WH-LRG	LRG	11-12 1/2	12-13 1/2
AKL-6334-BK-1XL	AKL-6334-WH-1XL	XL	13-14 1/2	14-15 1/2
AKL-6334-BK-2XL	AKL-6334-WH-2XL	2XL	15+	16+

Maximum Protection & Support

The Swede-O Strap Lok® ankle brace has external figure eight straps which allow you to easily control and customize the amount of support as well as simulate professional taping.

The exclusive dual eyelet lacing system is offset to create a "pulley effect", keeping the brace tighter longer. The top strap covers the laces and locks the figure eight straps in position. The full elastic back is comfortable against the Achilles tendon and helps reduce irritation and the potential for blistering.


The breathable tongue helps keep the foot cool and comfortable. The Strap Lok is backed by an industry leading 1 year warranty and constructed from a durable ballistic nylon which is comfortable thin, and lightweight.



Indications for Use

Use the Strap Lok after injury or to protect your ankle from future injury. This ankle brace helps prevent and relieve aches, pains, sprains, and other ankle conditions.

Available From:

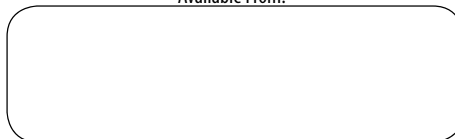

Specifications

Product Selection:

#AKL-6334 Swede-O Strap Lok®

Product Weight: .35 lb

Dimensions: 9" x 11.5" x 3"

Case Weight: 7 lbs

Case Dimensions: 11" x 14" x 13"

Case Quantity: 12

Materials: Polyester, nylon, polyurethane foam, neoprene, brass

