

M480


Embedding Cassettes

Made of acetal

Disposable plastic tissue cassettes are suitable for holding and identifying tissue samples in processing, embedding, and sectioning procedures. The cassettes fit securely in microtome chuck adapters. Cassette is totally resistant to the chemical action of histological solvents. These cassettes are designed to accept standard metal lids (cat.# M481) and will keep specimens in complete safety during processing. The slanted writing surface accepts markings easily, permitting sample identification throughout all stages of embedding and long afterwards when in archives. They are available in 11 colors. Each case contains 3 dispenser boxes of 500 cassettes.


Cat. #	Color	Qty/Cs	Cat. #	Color	Qty/Cs
M480-2	White	1500	M480-8	Tan	1500
M480-3	Pink	1500	M480-9	Gray	1500
M480-4	Green	1500	M480-10	Lilac	1500
M480-5	Yellow	1500	M480-11	Orange	1500
M480-6	Blue	1500	M480-12	Aqua	1500
M480-7	Peach	1500			

Cat. #	Description	Qty/pk
M481	Metal Process Cover	25

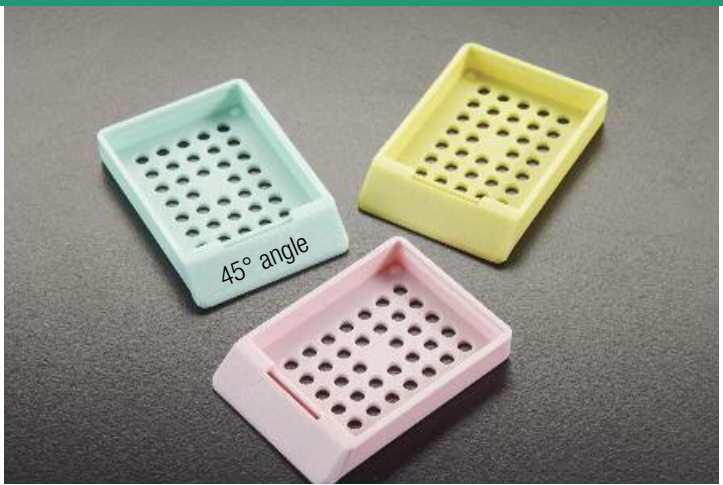

M481


Metal Lid for M480 Cassette

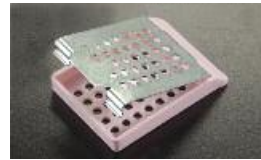

For IVD use



M480SL
For ThermoFisher
Printers

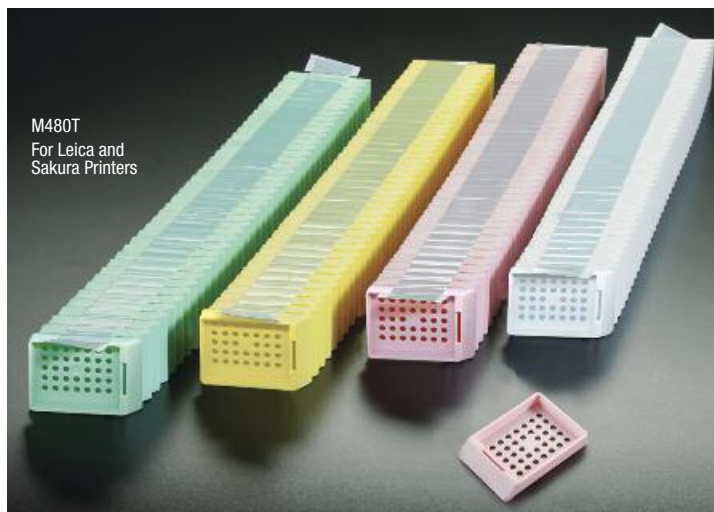

M480T
For Leica and
Sakura Printers

M480SL & M480T

Cassettes in QuickLoad™ Sleeves and Stacks

Made of acetal

The sleeved cassettes are especially made to be used with ThermoFisher printers. Cassettes with tape are to be used with Leica and Sakura Ink Jet printers. Molded from a special high density acetal polymer, they keep specimens safely submerged and are totally resistant to the chemical action of most solvents used in histology laboratories. The efficient flow-through slots maximize fluid exchange and ensure proper reagent drainage.

In Sleeves		In Stacks		Color
Cat. #	Qty/Cs	Cat. #	Qty/Cs	
M480-2SL	750	M480-2T	2000	White
M480-3SL	750	M480-3T	2000	Pink
M480-4SL	750	M480-4T	2000	Green
M480-5SL	750	M480-5T	2000	Yellow
M480-6SL	750	M480-6T	2000	Blue
M480-7SL	750	M480-7T	2000	Peach
M480-8SL	750	M480-8T	2000	Tan
M480-9SL	750	M480-9T	2000	Gray
M480-10SL	750	M480-10T	2000	Lilac
M480-11SL	750	M480-11T	2000	Orange
M480-12SL	750	M480-12T	2000	Aqua

Fluo Pink, Fluo Green and Fluo Yellow are also available. Minimum quantities apply. Contact us for more details.



Compatible with all
cassette printers

For IVD use