



M490 HISTOSETTE® I Tissue Processing / Embedding Cassettes

Made of acetal

Disposable plastic cassettes hold tissue specimens during the embedding process, as well as in a storage file. Molded from a special high density polymer, these cassettes keep specimens safely submerged in liquid and are totally resistant to the chemical action of histological solvents. **The efficient flow-through slots maximize fluid exchange and ensure proper drainage.** The one-piece integral lid eliminates the need for separate steel lids. Just snap apart and conveniently lock the lid into the base of the cassette. They can be opened or closed as often as necessary and they always relock securely without danger of specimen loss. Anterior writing area at a 30° angle. Not suitable for automated printers.

Each case contains 3 dispenser boxes of 500 cassettes.

For IVD use



M491 HISTOSETTE® I Biopsy Processing / Embedding Cassettes

Made of acetal

Disposable plastic cassettes similar to Series M490 but specially designed to hold biopsy specimens during the embedding process, as well as in a storage file. Anterior writing area at a 30° angle. Not suitable for automated printers.


Each case contains 3 dispenser boxes of 500 cassettes.

For IVD use

Tissue Cat. #	Biopsy Cat. #	Color	Qty/Cs
M490-2	M491-2	White	1500
M490-3	M491-3	Pink	1500
M490-4	M491-4	Green	1500
M490-5	M491-5	Yellow	1500
M490-6	M491-6	Blue	1500
M490-7	M491-7	Peach	1500
M490-8	M491-8	Tan	1500
M490-9	M491-9	Gray	1500
M490-10	M491-10	Lilac	1500
M490-11	M491-11	Orange	1500
M490-12	M491-12	Aqua	1500

Fluo Pink, Fluo Green and Fluo Yellow are also available. Minimum quantities apply. Contact us for more details.

Anatomy of a Histosette® I

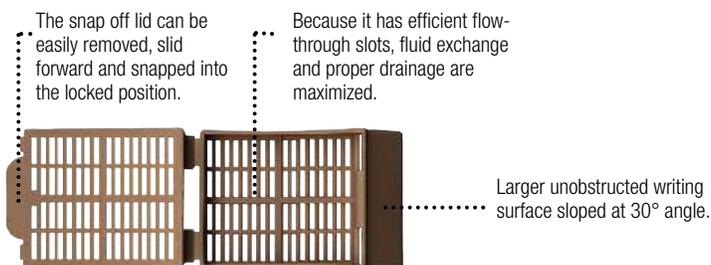

Because it has a back mounted locking device, it will never open during processing.

